



CHAIRSIDE®

Attachment Processing Material



#9566



Introducing — a new and unique Attachment Processing Material from ZEST Anchors, the manufacturer of the trusted LOCATOR® Attachment and LOCATOR Overdenture Implant (LODI) System. CHAIRSIDE® Attachment Processing Material was designed for ease of use and predictability when processing attachment components into overdentures, including ZEST's LOCATOR and SATURNO™ Denture Caps.

CHAIRSIDE Attachment Processing Material Key Attributes:

Cost Efficient

- No primer required for reduced chairtime
- Shorter mixing tip reduces material waste and cost per case*

Simple Application


- Ideal viscosity for accurate application
- Bonds to itself for efficient void fill and touch-ups
- Dual cure capability provides convenient self or light cure flexibility

Patient Friendly

- Low curing temperature ensures patient comfort
- Odorless and tasteless to increase patient acceptance

Interested in ordering?
Visit us on www.xpdent.com
or give us a call at
1.877.328.3965

Six Simple Steps For Chairside Processing LOCATOR® Attachments Into Overdentures

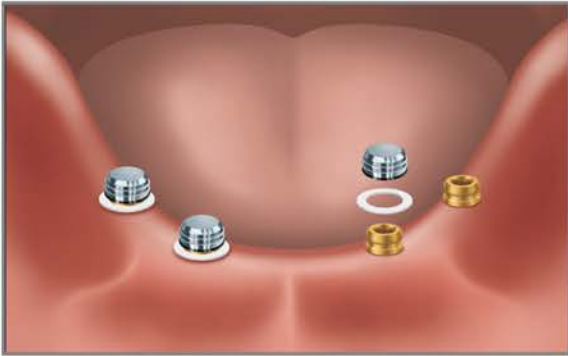


1. Place a White Block Out Spacer Ring around each abutment and press it down to the tissue. Place a Denture Cap with a Black Processing Male onto each abutment and press down firmly to ensure engagement of the abutment.



2. Apply fit check marking paste inside the denture. Insert it into the mouth over the Denture Caps. Mark the areas where the denture will need to be relieved to allow space for the caps to be picked up.



3. Relieve the areas marked with an acrylic bur. Try in the denture to verify that the Denture Caps are not in contact with the acrylic in any area. Cut lingual/palatal vent windows in the denture to visualize full seating and for excess material to vent. Cut an undercut below the intaglio surface of each Denture Cap relief area for mechanical retention.

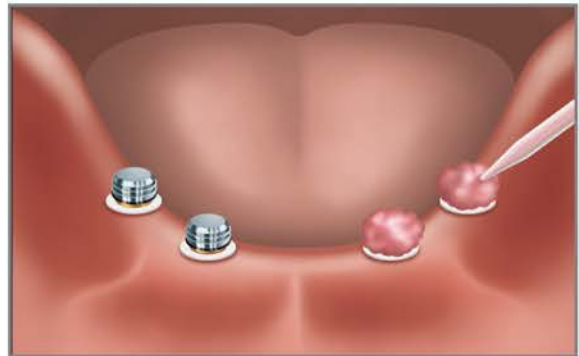

4. Dry the Denture Caps. Apply a small amount of CHAIRSIDE Material around the circumference of each cap. Place CHAIRSIDE Material into the relief areas of the denture and seat it over the caps and onto the tissue. Have the patient close into light occlusion and hold while the CHAIRSIDE Material sets. The set time of the material from the time the material is mixed to the time it is set is five (5) minutes. The working time for the material is two (2) minutes, fifteen (15) seconds.



5. Disengage the denture from the abutments and remove from the mouth. Verify that the Denture Caps have been securely processed into the denture. Fill any voids, light cure and polish the denture. The material will bond to itself and will cure within 30 seconds with light application.



6. Place the denture in the mouth and press down to engage the males on the attachment system.