

**WE'VE MADE
FIXED FULL-ARCH
RESTORATIONS
A SNAP.**

**NO SCREWS. NO CEMENT.
NO COMPROMISES.**

What if screws or cement were
no longer needed to attach a fixed
prosthesis to the abutments?

What if you could shorten and
simplify the treatment appointment?


What if the patient's esthetics and
comfort were maximized?

These questions were at the foundation of the design and development for the unique and innovative **LOCATOR F-Tx[®] Fixed Attachment System** from Zest Dental Solutions. By leveraging our many years of expertise with the LOCATOR[®] Attachment System and our understanding of full-arch solutions for edentulous patients, we strived to create a simpler, more efficient system for **FIXED** full-arch implant restorations.

Fixed for the Patient. Easily Removed by the Clinician.

Welcome to LOCATOR F-Tx
- No Screws. No Cement. No Compromises.

INTRODUCING THE FIRST FIXED FULL-ARCH REST

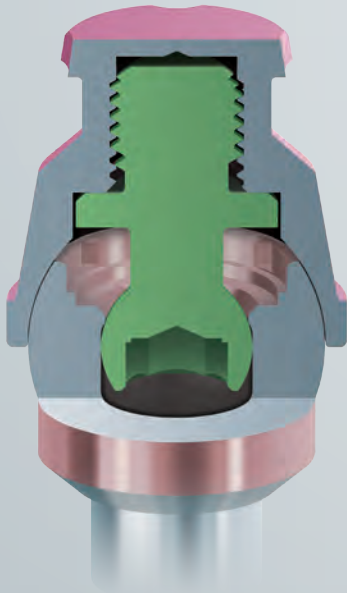

SECURE FIT DESIGN


This revolutionary, patent-pending retention system works similar to a ball and socket, allowing the Denture Attachment Housing to securely seat into place and then pivot to the desired position. Once in place, it's fixed for the patient, and can easily be removed by the clinician during hygiene and maintenance visits.

EASILY AND QUICKLY SEAT AND REMOVE THE PROSTHESIS


One tool for both seating and removal of the LOCATOR F-Tx prosthesis.



- Simple and confident seating of the LOCATOR F-Tx Prosthesis.

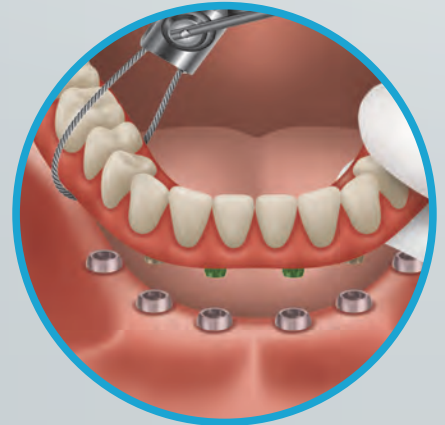

- Patient friendly and quick removal of the LOCATOR F-Tx Prosthesis.

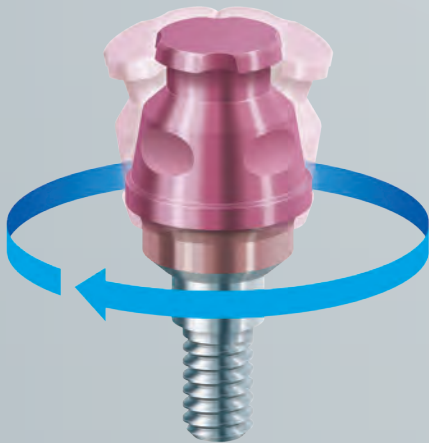

Denture Attachment Housing pivots 20 degrees in any direction allowing for the Housing to be positioned into the ideal location for the prosthesis.

SIMPLIFIED ANGLE CORRECTION/STRESS FREE PASSIVE FIT

- The LOCATOR F-Tx Abutment features a unique, spherical coronal geometry which allows the Denture Attachment Housing to rotate in any direction and correct up to 20 degrees from a common vertical eliminating the need for angled abutments.
- Chairside processing procedures at final prosthesis delivery ensures a stress free passive fit.

LOCATOR F-Tx®.

PROVISIONALIZATION WITHOUT COMPROMISE.




DENTURE ATTACHMENT HOUSING

- Denture Attachment Housing is threaded internally to accept a PEEK Retention Ball that snaps into the LOCATOR F-Tx Abutment.
- Features an anodized pink finish for improved esthetics.
- Aggressive grooves and flats limit vertical and rotational movement.
- Denture Attachment Housing is passively picked-up in the prosthesis via a chair side technique.



RETENTION/PROCESSING BALLS

- PEEK Retention Balls are available in low, medium and high retention levels based on the needs of the case.
- A Processing Ball comes pre-inserted with the Denture Attachment Housing, and an additional Processing Ball is included, both are used for provisionalization and laboratory procedures.

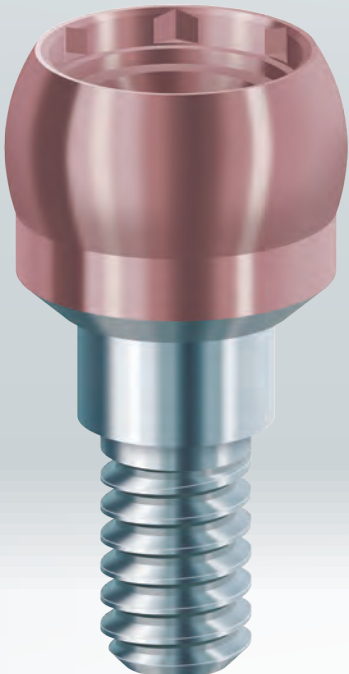

LOW



MEDIUM



HIGH



PROCESSING

ABUTMENT

- Unique, spherical coronal geometry allows the Denture Attachment Housing to correct up to 20 degrees from a common vertical and pivot in any direction to position the Housings into the ideal location for the prosthesis.
- DuraTec® Coating provides a hard, smooth and wear-resistant abutment exterior with an esthetically pleasing gingival tone.

ALL-IN-ONE PACKAGING


- LOCATOR F-Tx features all-in-one packaging that is sterile and includes everything you need. Abutment (with cap to deliver the abutment to the implant site), Denture Attachment Housing with pre-inserted Processing Ball, an extra Processing Ball, as well as one Blue (**Low**), Tan (**Medium**), and Green (**High**) Retention Ball.



LOCATOR F-Tx[®]

A NEW WAY TO THINK ABOUT FIXED FULL-ARCH RESTORATIONS.

LOCATOR F-Tx IS A SIMPLIFIED, TIME SAVING SOLUTION for fixed full-arch restorations with no compromise to prosthesis strength or esthetics. Optimized for efficiency and chair time savings compared to conventional screw-retained systems, it features a novel, patent-pending attachment that eliminates the need for sub-gingival cement or screw access channels. An exceptional solution for immediate provisionalization, the fully integrated system accommodates the full range of final restorative options including all acrylic, metal-reinforced acrylic, PFM, and all zirconia frameworks. **Truly unique and cutting-edge, LOCATOR F-Tx is sure to deliver satisfaction for both the patient and clinician.**

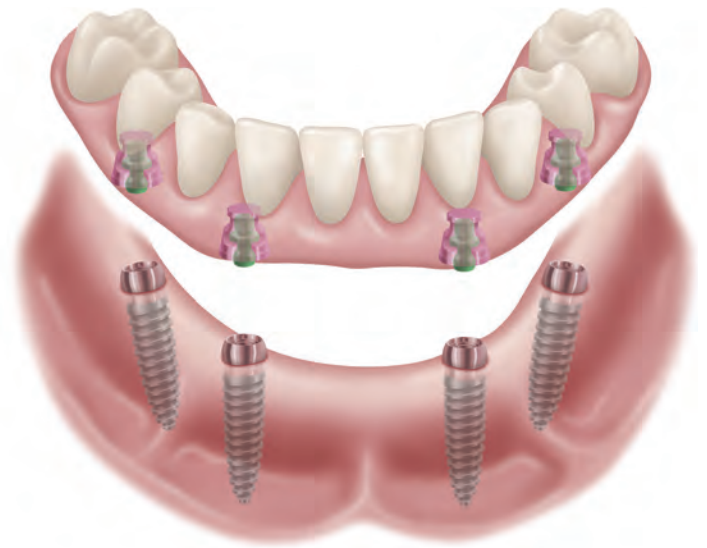

SIGNIFICANT TIME SAVINGS

Significant chair time savings over conventional techniques for fabrication of an immediate provisional prosthesis. Patient visits are also expedited and simplified with quick, easy removal of the prosthesis during follow-up visits for final impressions and try-ins, as well as future recall and maintenance visits.

SIMPLICITY

Easily seated attachments and minimal components streamline procedures, while the unique design of the LOCATOR F-Tx Fixed Attachment System accommodates scenarios up to 20 degrees from a common vertical eliminating the need for angled abutments.

STRENGTH AND ESTHETICS

No cylinders, screw access channels, and screw heads to compromise esthetics or weaken the prosthesis. The ability of each Denture Attachment Housing to pivot on the Abutment up to 20 degrees in any direction allows the housing to be positioned into the ideal location for the prosthesis.



**SCREW RETAINED PROSTHESIS
WITH SCREW ACCESS CHANNELS**



**LOCATOR F-Tx PROSTHESIS
MAINTAINS INTEGRITY AND ESTHETICS**



To learn more about the simplicity, time savings, and unique attachment design of the new LOCATOR F-Tx® Fixed Attachment System please visit us at
www.zestdent.com/ftx



2875 Loker Avenue East
Carlsbad, CA 92010 USA
www.zestdent.com

Toll Free: 1.800.262.2310
Direct: 1.442.244.4835
Fax: 760.743.7975
Email: sales@zestdent.com

DISTRIBUTED BY: