



LOCATOR® ATTACHMENT SYSTEM

The worldwide standard of care for overdentures!




The **SELF-ALIGNING** feature of the LOCATOR Attachment allows a patient to easily seat their overdenture without the need for accurate alignment of the attachment components.



NEW!
8 lb Retention
Level Available

Interested
in ordering?
Visit us on
www.xpdent.com
or give us a call at
1.877.328.3965

Locator Retention Levels

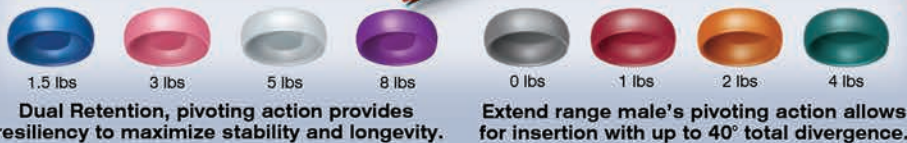

Implant Attachment



Bar Attachment



Root Attachment

ZEST
ANCHORS

Your East Coast Supplier for Zest Locator Attachment System


XPdent Corporation Sales: **1.877.328.3965** Tech Support: **1.877.329.3965** Web: www.xpdent.com

Locator® Implant Attachment Quick Reference Guide

The LOCATOR Implant Attachment is available in stock to fit all major implant systems. LOCATOR Implant Abutments for overdenture applications are available in numerous tissue cuff heights. Choose the corresponding abutment tissue cuff height that exactly equals the tissue measurement, or the next higher size. With proper tissue cuff height selection, the 1.5mm working portion of the attachment will remain above the gingival tissue.



LOCATOR
Implant Abutment and titanium Denture
Cap with clear nylon Replacement Male

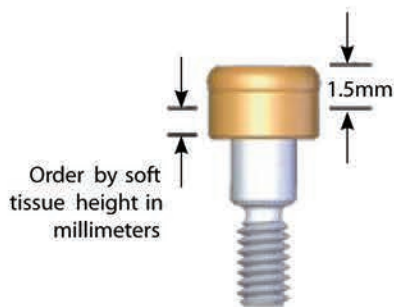

Order by soft
tissue height in
millimeters



LOCATOR Implant Abutments
Many Tissue Cuff Heights are available.

The LOCATOR Core Tool (8393) is necessary for abutment placement and nylon male placement & removal; includes Male Removal Tool, Male Seating Tool, and Abutment Driver, shown with optional Abutment Holder Sleeve (8394). Various torque wrench connection types of insert drivers are available to achieve optimum 30N-cm of torque.



LOCATOR Core Tool
with optional Abutment Holder Sleeve

LOCATOR Insert Driver
(square drive connection)



LOCATOR Insert Driver
(latch-type connection)



LOCATOR
Male Processing Package



LOCATOR Extended Range
Male Processing Package

LOCATOR Male Processing Packages are available for parallel implants (8519-2, 8519-10) and divergent implants (8540-2, 8540-10). A full range of retentive males are included with each denture cap to allow personalized retention for each specific patient (nylon male retention selection depends upon patient preferences and number of implants placed into function).

LOCATOR Male Processing Package (8519) includes Denture Cap with Black Processing Male, Dual Retentive Replacement Males: Clear, Pink, Blue, and Block-Out Spacer. Offered in packs of two (2) or packs of ten (10).

LOCATOR Extended Range Male Processing Package (8540) includes Denture Cap with Black Processing Male, Extended Range Replacement Males: Green, Orange, Red, and Block-Out Spacer. Offered in packs of two (2) or packs of ten (10).



Dual Retention 0-10 degrees
(8559) Purple: 8.0 Lbs (3629 Grams)
(8524) Clear: 5.0 Lbs (2268 Grams)
(8527) Pink: 3.0 Lbs (1361 Grams)
(8529) Blue: 1.5 Lbs (680 Grams)



Extended Range
Angle Retention 10-20 degrees
(8547) Green: 3-4 Lbs (1361-1814 Grams)
(8915) Orange: 2.0 Lbs (907 Grams)
(8548) Red: 0.5 - 1.5 Lbs (226-680 Grams)
(8558) Gray: 0 (zero) Retention

Parallel posts (8517) and the angle measure guide (9530) are useful to determine the appropriate nylon male selection based upon angulation.



LOCATOR Parallel Post



Angle Measurement Guide